

Infinity 1046

窒素発生装置

品番: 3302084




説明

最大20 L / minの窒素と70 L / minの乾燥空気を生成するInfinity 1046は、シングル/デュアルソース構成にかかわらず、PerkinElmer QSight Triple Quadをサポートするように特別に設計されています。Infinityシリーズは、コンプレッサーを使用しないガス発生器であり、効果的に動作するためにクリーンな室内空気源を必要とします。可動部品が非常に少ないため、これらの発電機は効果的に静音で動作し、継続的な高性能と製品寿命の延長を確実にするために最小限の注意のみが必要です。



Generator Details

- ガスの種類: 窒素, 空気
- 統合型コンプレッサ: 番号

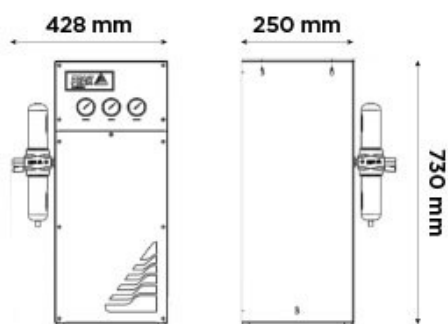

主な機能

- 最大20 L / minの窒素（乾燥ガス）と最大70 L / minの空気（ソースガス）を供給
- シングルソースとデュアルソースの両方のQSightモデルを提供可能
- 必要に応じて、最適なパフォーマンスで24時間365日稼働
- ライフタイムのランニングコストが低く、非常に経済的な窒素と空気の供給源
- 最小限の設定が必要
- 完全に無音の動作
- 12か月の包括的なオンサイト保証

技術仕様

	Infinity 1046 窒素発生装置
ガスの種類	窒素
最大ガス流量:	20L/min
最大出口圧力:	110psi/7.6bar
ガスの種類	空気
最大ガス流量:	70L/min
最大出口圧力:	110psi/7.6bar
最大出口圧力:	110psi/7.6bar
ガス出口の設定:	2 x 1/4" BSPP
最低入口圧力:	120psi/8.3bar
最高入口圧力:	145psi/10bar
Min Air Inlet Flow:	180
最小入口空気品質:	ISO8573 - 1:2010 Class 1.4.1
起動時間:	60 mins
消費電力:	0
電圧:	n/a
周波数:	n/a
電流:	n/a
熱出力:	n/a
動作時の最高温度:	30°C / 86°F
分子:	<0.01µm
大きさ (高さx幅x奥行き) ミリメートル	730 x 422 x 250 mm
大きさ (高さx幅x奥行き) インチ	28.7 x 16.6 x 9.8 inches
装置の重量	23kg / 50.6lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

