

CG15L Zero Air

校正ガスの連続供給

品番: 43-1010 (230V) - UK | 43-1020
(230V) - EU | 43-1030 (110V) - US




説明

大気放出アプリケーションなどの多くの実験室アプリケーションでは、校正ガスが必要です。これらの要件を満たすために、Peak Scientificは、ゼロエアモデルよりも高いパフォーマンスを必要とするアプリケーション向けに、二酸化炭素、メタン、および水分を空気から除去することによりゼロエアを生成するCG15Lを製造しました。



Generator Details

- ガスの種類: 空気
- 統合型コンプレッサ: 番号

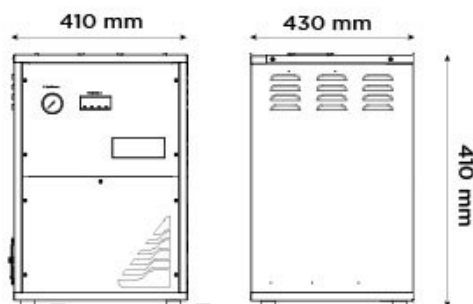

主な機能

- クリーンでドライな実験室グレードの空気を生成し、ガス発生器と機器への最適な供給を保証します
- 校正ガスの継続的な供給
- 連邦レジスター仕様40CFR 86.114-79の要件に適合
- 可動部品が少ないため、メンテナンスがほとんど不要で、発電機の寿命が長くなります
- 完全に無音の動作
- Precision Air Compressorでサポート可能
- 最小限の設定で12か月の包括的なオンサイト保証が必要

技術仕様

	CG15L Zero Air
ガスの種類	空気
最大ガス流量:	15L/min
最高純度:	<0.1ppm
ガス出口の設定:	1x 1/4" BSP Female
最低入口圧力:	100psi/6.9bar
最高入口圧力:	125psi/8.6bar
起動時間:	30 mins
消費電力:	1100/540 W
電圧:	120 / 220 V ±10%
周波数:	50 / 60 Hz
電流:	2.6 / 5.2 A
認証:	CE
大きさ (高さx幅x奥行き) ミリメートル	410 x 410 x 430 mm
大きさ (高さx幅x奥行き) インチ	16.1 x 16.1 x 16.9 inches
装置の重量	58kg / 127.6lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

